




Providing Relief for Patients with Dry Eye or Lacrimal Disorders

Parasol® | Extend® Punctal Occlusion



**Permanent
Parasol®**

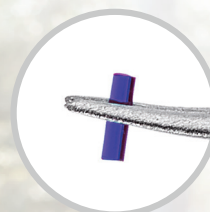

**Long-Term
EXTEND 180®**



**Mid-Term
EXTEND®**



**Partial
Micro Flow™**



**Short-Term
Collagen**

Permanent Occlusion

Permanent occlusion is indicated for chronic dry eye conditions. All occluder products are made of medical grade silicone and are reversible.

Insertion

The hollow nose collapses on insertion, eliminating the need for punctal dilation in most cases. It's straight-forward for you to insert and more comfortable for your patients.

Straight Forward Sizing

Two sizes, small and medium, fit virtually all patients. X-Small and large sizes are available for customization.

Retention

Parasol is designed for excellent retention, as demonstrated by an independent clinical study.¹

Parasol® Size Comparison Chart



Traditional Plugs

0.2 mm	0.3 mm	0.4 mm	0.5 mm	0.6 mm	0.7 mm	0.8 mm	0.9 mm +
--------	--------	--------	--------	--------	--------	--------	----------

Parasol Plugs

X-SMALL	SMALL	MEDIUM	LARGE
---------	-------	--------	-------



Parasol® Pre-Loaded Punctal Occluders Sterile

Occluders are packaged in a convenient 1 pair tray, are sterile and pre-loaded on a single-use inserter.



A14-201



QTY: 2/Bx

Parasol® Pre-Loaded Punctal Occluders Sterile	SKU #
X-Small, 0.25 - 0.35 mm	A14-200
Small, 0.35 - 0.65 mm	A14-201
Medium, 0.6 - 0.85 mm	A14-202
Large, 0.9 mm+	A14-203

Parasol® Quickload Punctal Occluders Sterile

Occluders are sterile and pre-mounted on a stand for Quickload convenience. Reusable inserters are sold separately.

QTY: 16/Bx



Parasol® Quickload Punctal Occluders Sterile	SKU #
Small, 0.35 - 0.65 mm	A14-201QL
Medium, 0.6 - 0.85 mm	A14-202QL
Large, 0.9 mm+	A14-203QL
Combo small and medium	A14-205QL

¹ C. McCabe, MD, PhD, FACS Plugs Reduce Dry-Eye Symptoms, Improve Vision. November 2009, Review of Ophthalmology.

Parasol® Bulk Punctal Occluders Non-Sterile

Occluders are non-sterile and packaged 2 occluders per poly bag. Reusable inserters are sold separately.

QTY: 20/Bx



Parasol® Bulk Punctal Occluders Non-Sterile	SKU #
Small, 0.35 - 0.65 mm	A14-201BN
Medium, 0.6 - 0.85 mm	A14-202BN
Large, 0.9 mm+	A14-203BN
Combo small and medium	A14-209BN

Plus Pre-Loaded Punctal Occluders Sterile

The Plus design is the result of a micro-molding manufacturing process using a soft medical grade silicone.

- A sharply angled nose for insertion
- An ultra-smooth surface finish
- Delicate contoured edges
- Designed for consistent retention

Occluders are packaged in a convenient 1 pair tray. Each occluder is sterile and pre-loaded on a high quality disposable inserter.

QTY: 2/Bx



A14-204

Plus Pre-Loaded Punctal Occluders Sterile	SKU #
Small, 0.4 - 0.55 mm	A14-204
Medium, 0.55 - 0.7 mm	A14-205
Large, 0.7 - 0.85 mm	A14-206
X-Large, 0.85 mm+	A14-207

Plus Punctal Occluders Non-Sterile Bulk

Occluders are supplied non-sterile and have 2 occluders per poly bag, with 10 bags/pairs per box, unless otherwise noted. Reusable inserters are sold separately.



Plus Punctal Occluders Non-Sterile Bulk	SKU #
Small, 0.4 - 0.55 mm	A14-204BN
Medium, 0.55 - 0.7 mm	A14-205BN
Large, 0.7 - 0.85 mm	A14-206BN
X-Large, 0.85 mm+	A14-207BN
Large, 0.9 mm+	A20-203BN

Mid-Term Occlusion

Extend® and Extend® 180 Synthetic Punctal Occluders

Provides occlusion for up to 3 months with Extend® or 6 months with Extend® 180, for patients experiencing short-term dry eye symptoms. These are ideal for the transient dry eye associated with LASIK or other ocular surgeries. Can also be used for allergy sufferers, contact lens discomfort, climate change or to sustain drug contact time and effectiveness.

QTY: 20/Bx



A13-204



A13-2044

Extend® Punctal Occluders	SKU #
Variety Pack	A13-200
0.2 mm x 2 mm	A13-202
0.3 mm x 2 mm	A13-203
0.4 mm x 2 mm	A13-204
0.5 mm x 2 mm	A13-205

QTY: 20/Bx

Extend® 180 Synthetic Punctal Occluders	SKU #
Variety Pack	A13-2000
0.3 mm x 2 mm	A13-2033
0.3 mm x 2 mm, Single Pair	A13-2033SP
0.4 mm x 2 mm	A13-2044
0.4 mm x 2 mm, Single Pair	A13-2044SP
0.5 mm x 2 mm	A13-2055
0.5 mm x 2 mm, Single Pair	A13-2055SP

Partial Occlusion

Micro Flow™ Punctal Occluders

Designed to provide a straight-forward, reversible means of partial occlusion. Micro Flow™ Punctal Occluders support lacrimal drainage care. Suggest to use when total occlusion is not ideal.

QTY: 2/Bx



A14-401

Micro Flow™ Punctal Occluders	SKU #
Sterile, Small, 0.4 - 0.55 mm	A14-401
Sterile, Medium, 0.55 - 0.7 mm	A14-402
Sterile, Large, 0.7 - 0.85 mm	A14-403

QTY: 10/Bx

Micro Flow™ Punctal Occluders	SKU #
Non-Sterile, Bulk, Small, 0.4 - 0.55 mm	A20-401BN
Non-Sterile, Bulk, Medium, 0.55 - 0.7 mm	A20-402BN
Non-Sterile, Bulk, Large, 0.7 - 0.85 mm	A20-403BN

Temporary Occlusion

Collagen Synthetic Punctal Occluders

These temporary occluders provide occlusion for 7-10 days. They are ideal for diagnostic use for determining if permanent punctal occlusion is needed for your patient. 6 occluders in each foam sleeve, with 10 foam sleeves per box.

Qty: 60/Bx



A12-103

Collagen Synthetic Punctal Occluders	SKU #
0.2 mm x 2 mm	A12-102
0.3 mm x 2 mm	A12-103
0.4 mm x 2 mm	A12-104
Variety Pack	A12-200

Single-Use Punctal Dilators and Probes

Punctum Dilator, Fine/Medium

0.4 mm x 21.5 mm Tip,
0.6 mm x 18 mm Tip,
77 mm Handle



Qty: 20/Bx
Brand: Malosa®

Description	SKU #
Punctum Dilator, Fine/Medium	MMSU1181S

Punctum Dilator/Seeker

0.45 mm x 21.5 mm Tip,
5 mm Pointed tip,
18 mm Straight shaft, 77 mm Handle




Qty: 20/Bx
Brand: Malosa®

Description	SKU #
Punctum Dilator/Seeker	MMSU1193S

Visitec® Punctal Dilator and Probe, Khan

Dual function instrument allows for single handed transition from dilating to probing without release of lateral lid traction.

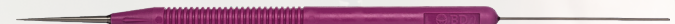

Qty: 3/Bx
Brand: Visitec®

Description	SKU #
25 mm/1 in medium taper, Blunt, Wilder style dilator mated with a size 0 Bowman lacrimal probe, 0.7 mm x 45 mm/22G x 1 3/4 in	585124

Other Related Occlusion Options

Curved Tying Forceps

0.25 mm Tips, 5 mm Curved tying platforms, 80 mm Handle



Qty: 20/Bx
Brand: Malosa®

Description	SKU #	SKU # Clips
Curved Tying Forceps	MMSU1210S	MMSU1210CS

Lacrimal Cannula



Qty: 10/Bx
Brand: Visitec®


Description	SKU #
● 25G x 1/2 in/0.5 x 13 mm	585068

Lacrimal Cannula



Qty: 10/Bx
Brand: Visitec®

Description	SKU #
● 26G x 1 1/8 in/0.45 x 28 mm	581276

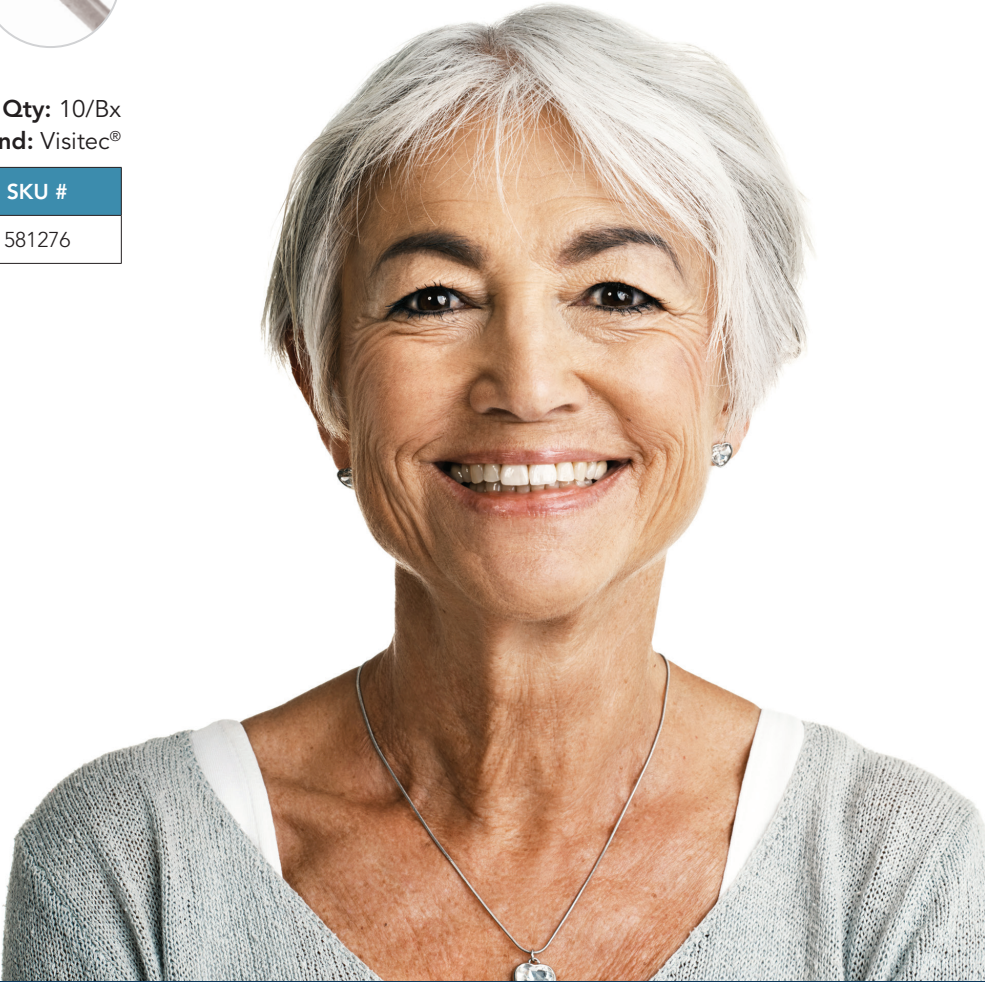

FOR COMPLETE PRODUCT INSTRUCTIONS,
PLEASE REFER TO THE BVI IFU OR PRODUCT INSERT

Tel: 866-906-8080
Email: customersupport@bvimedical.com

BVI and all other trademarks (unless noted otherwise) are property of BVI. BVI ©2021